



STORAGE

Instruction Manual

Controlled Medical Device: Electrical Massager for Home Use

HAND CARE AX-HXL180




※Design specifications may have been partially changed. Thank you for your understanding

Thank you very much for purchasing AC-HXL180. Before use, read this instruction manual and use the product accordingly. Please keep this manual with care after reading it.

ATEX


The names and functions of the parts

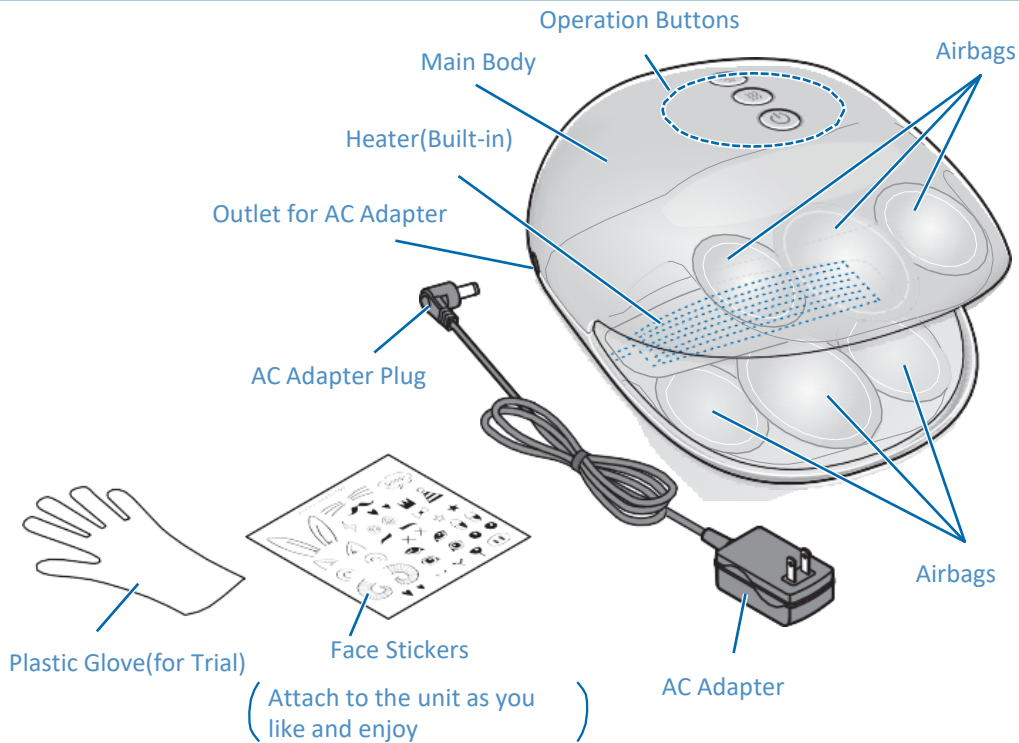

How to Use

CAUTION Precaution before Use


- Check whether there are any sewing tears on the massaging area. Do not use if it is torn.
- Make sure the switch and timer are functioning properly before use.
- After long periods time of disuse, reread the instruction manual well and make sure that the product functions normally and safely prior to use.

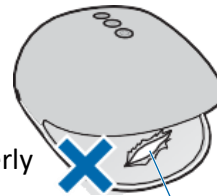

Massaging Area
Please check whether there are any sewing tears.


1 Connecting AC adapter and body.

Please insert them straight into the root securely.

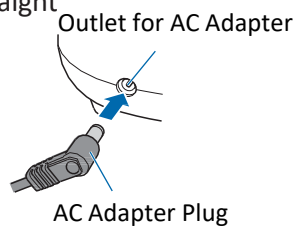


2 Plug in the AC adapter to the outlet

Make sure the adapter is plugged in correctly and securely.

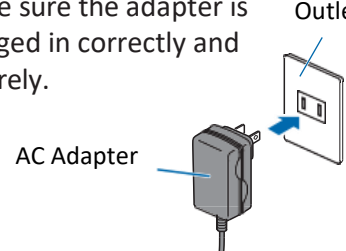

3 Insert your hand in the device.

Open up the hand you wish to use and insert four fingers excluding the thumb into the finger holes of the device.

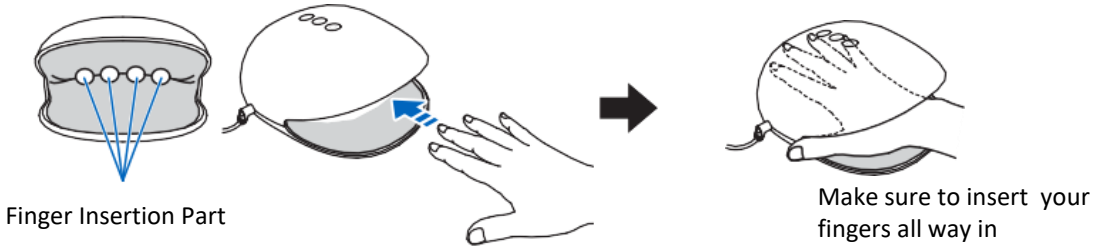


Example of Use




CAUTION

Always turn the power off when inserting and removing hands from device.

Left Hand

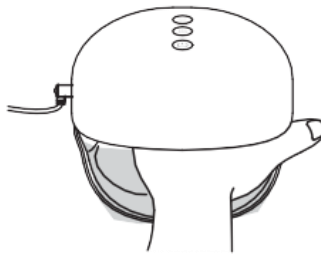

Right Hand

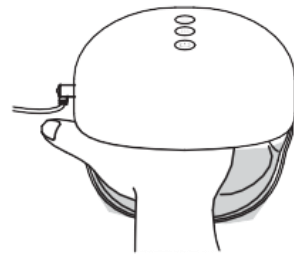


To massage the thumb, turn your hand sideways and insert only the thumb into the device.

(Example of Using Thumb Massage)

Left Hand




Right Hand

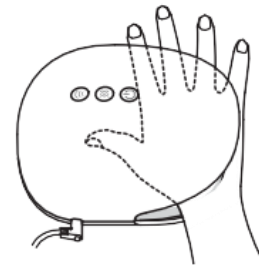

(How to Use Plastic Glove)

If you wish to apply hand cream when using plastic gloves, use what you normally use.

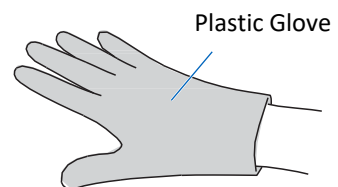

The plastic gloves that come with the product are disposable and to be used as a sample. Once you've used up the plastic gloves that come with the product, use ordinary plastic gloves sold in stores. Do not reuse for sanitary reasons.

If your hand is pushed out of the device while wearing a plastic glove, hold the device down lightly with the other hand.


4


Push the power button to start the massage.

(Start massage)


Press  the power button to start the full massage course.


(Choose massage course)

By pressing  the power button the unit mode changes between the Full Course, Fingertip Course, and Power Off.

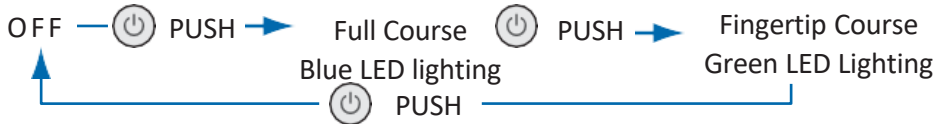


(Choose the strength of massage)


By pressing the strength button , the unit mode changes between strong and weak.




At the start, it is set to weak. By pressing the Strength button , the unit mode changes between strong ↔ weak.

(Start the heater)

By pressing the heater button , the unit Mode changes between on ↔ off.

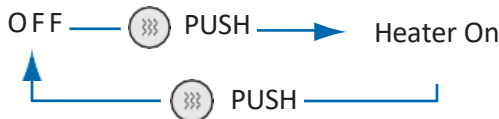At the start, it is set to off.




-  Strength Button
Weak: no light
Strong: Lighting
At the start, it is set to weak
-  Heater Button
✕red light when starting
-  Full course
✕blue light
Fingertips course
✕green light

- Air massage and heater can be used separately or at the same time.
- When about 10 minutes elapse from the time either the air massage or heater have started, the timer feature will activate and turn them off automatically.

Maintenance

- Please wipe the product with dry, soft cloth.
- For stains on the main body, dip a piece of soft cloth in a pH neutral detergent diluted by lukewarm water, wring it dry and wipe the stains off, then wipe the product with a piece of dry cloth.
- In order to prevent discoloration, do not use thinners or benzene.
- Wipe the metal part of the power plug once a year to prevent the insulation deterioration.

 WARNING	Be sure to turn off the power before disconnecting the AC adapter from the outlet.
---	--