



製造番号

HANDY MASSAGER *Rirabbit*

Instruction Manual

For storage

Medical device: home electric massager


Lourdes Handy Massager "Rabbit"
AX-KXL3400



*Please understand some designs and specifications may have changed.

Thank you very much for purchasing AX-KXL3400. Before use, read this instruction manual and use the product accordingly. Please keep this manual with care after reading it.

Name of the parts

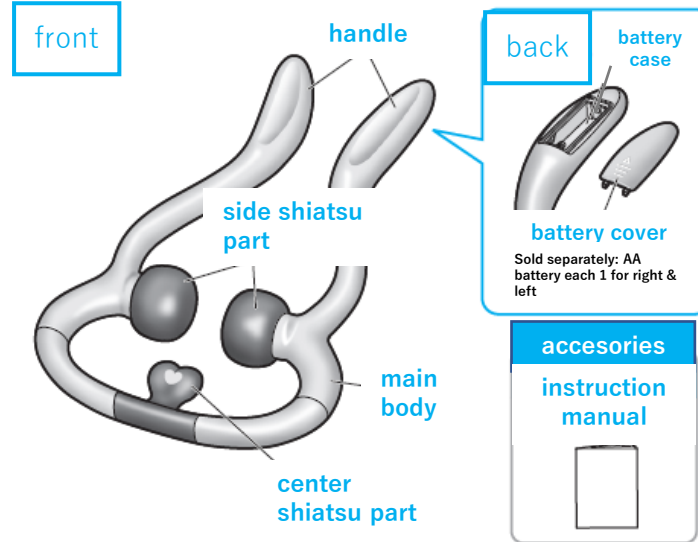


Preparation before use

How to put batteries

STEP 1 Remove battery cover.



STEP 2 Insert batteries according to + / - of battery case.

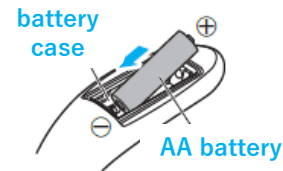

STEP 3 Attach battery cover.



How to use



Prohibited


- *Do not soak, use in bathroom or wash with water.
- *Do not use on parts other than the neck.
- < It may cause an accident or injury.>
- *Do not open the left & right handles with too much force.
- < It may cause an accident or failure.>

Compatible with waterproof function, IPX5.

*IPX5: There is no effect on water jets from all directions.



As figure, hold the handle with both hands, face up.

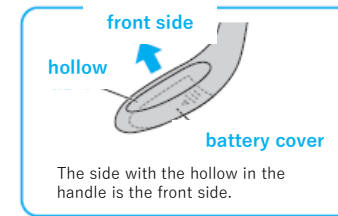

Move the handle in the direction of the arrow to massage.

Please massage with weak stimulation at the beginning of use.



The massage ends when the side shiatsu parts are released from the body.



CAUTION

The hair may be pinched by the shiatsu part.